



servvox

voice restoration and amplification solutions



KAPITEX

The Specialists in Laryngectomy,
Tracheostomy, Dysphagia and
Airway Management

Servox® digital

ELECTRONIC LARYNX



The latest model from the leading electronic larynx Company in the world yet again sets new standards of excellence.

Sophisticated electronics, to give maximum quality of sound and rehabilitation over many years of use - the hallmark of Servox.

features

- Individually programmable buttons - enabling the possibility of programming both buttons with different volume and frequency settings.
- In addition to the basic manual settings, the digital electronics make it possible to connect the speech aid to your PC via a serial interface. Special software enables you to make and save a multitude of individual settings. A useful tool for Speech Therapists, when using the Digital Data Communication Kit.
- A flashing red light on the Servox digital indicates that the battery is on reserve and requires charging. An elaborate energy management system protects the battery from being completely discharged while speech duration is increased by reduced energy consumption.
- The intelligent charging technique offers the simplest, safest and quickest way to charge your reliable Servox batteries. The battery can be kept in the Servox digital during charging. A replacement battery can also be charged at the same time.

KAPITEX HEALTHCARE - THE SPECIALISTS IN VOICE RESTORATION & AMPLIFICATION

Kapitex Healthcare brings you a superb choice of electronic larynxes and amplifiers. It allows the speech therapist to prescribe the unit best suited to the individual need.


Our new range has been carefully chosen by the specialists at Kapitex to enhance the voice rehabilitation of the laryngectomee and those with weak or whisper voice...enabling speech for everyday living.



All larynxes supplied by Kapitex are guaranteed for 24 months

A new dimension in communication offering a near full range of sounds

- The Servox digital has many user-friendly details, one of these being that the user can change or remove the carrying cord. If the carrying cord has been removed, the opening can be closed using the cord protector.
- The dark blue nylon hip bag gives the Servox digital a compact and up-to-date look. A separate compartment makes it easy to safely carry a replacement battery. The Velcro closure on the rear side allows you to adjust and attach the bag to suit different belt widths.
- Using the latest technology, perfectly constructed workmanship and a high-quality titanium housing, we can guarantee long durability, safety and optimal speech readiness at all times.
- The high value of the device is emphasised by the classic design and attractive colour scheme.

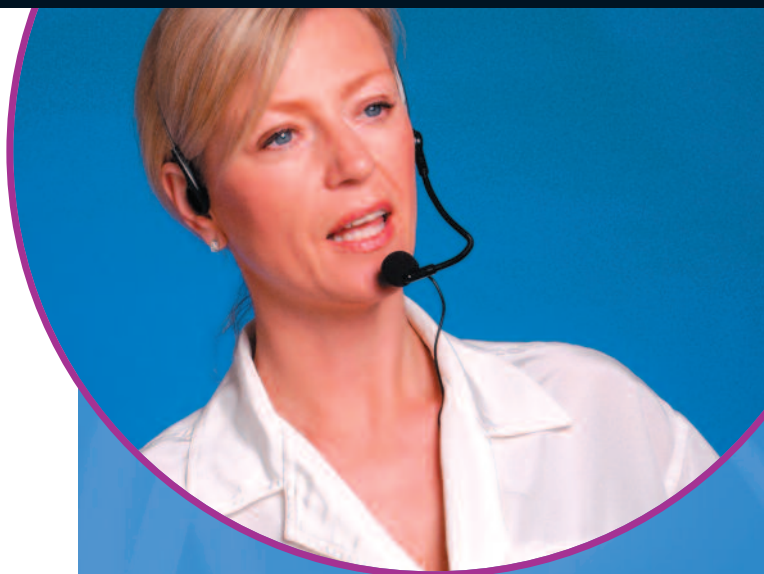


Digital Electronics with Serial Interface

In addition to the basic manual settings, the digital electronics make it possible to connect the speech aid to your PC via a serial interface. Special software enables you to make a multitude of individual settings.




Processor-Controlled Charger

This intelligent charging technique offers the simplest, safest and quickest way to charge batteries. The battery can be kept in the Servox digital during charging. A replacement battery can also be charged at the same time.

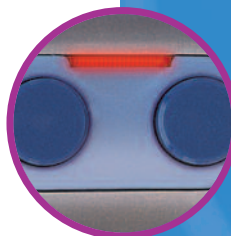

Optical Battery Control

A flashing red light on the Servox digital indicates that the battery is on reserve and requires charging. An elaborate energy management system protects the battery from being completely discharged while the speech duration is increased by reduced energy consumption.

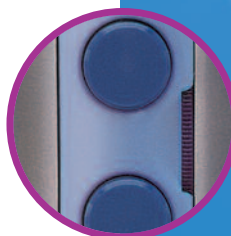

Individually Programmable Buttons

The Servox digital is the only speech aid which offers the possibility of programming both buttons with different volume and frequency settings. This makes it possible to pre-set the device for use in certain situations. The software also offers additional options for adjusting the tone to match the natural speech melody.

Technical Specification

Supplied with:

- 2 NiMH 8.4V batteries
- Battery charger
- Carrying pouch
- User instructions

Length: 11.7cm
Weight : 150gms

Manufactured in
Germany

Servox[®] PA5 POWER

PORTABLE VOICE AMPLIFIER

PA5 POWER

PA5 power is simply the finest portable voice amplifier. It can boost your volume by as much as 18 decibels. It consists of a rechargeable 'waist pack' amplifier and speaker unit along with an extremely comfortable headset microphone. You simply charge the PA5 power overnight with the included charger.

features

- Lightweight Portable Voice Amplifier (1lb 6oz including batteries)
- Built in rechargeable, environmentally friendly Nickel metal-Hydride (NiMH) batteries provide up to 14 hours of use. Recharging is easy and can be done overnight or as quickly as five hour charge. (For added insurance, you can even substitute AA batteries if you forget to charge the PA5 power.)
- Ergonomic design for comfortable use all day long
- Simple to use
- Adjustable volume
- Super portable



Technical Specification

Supplied with:

- Rechargeable batteries
 - Battery charger (Order No. LACSA 0205)
 - Headset microphone
- Manufactured in Germany

Length: 200mm
Height: 70mm
Width: 60mm
Weight : 800g
Power Output: 5w

www.kapitex.com

DEPEND ON KAPITEX QUALITY AND SERVICE

Kapitex Healthcare provide a wide range of products dedicated to laryngectomy and tracheostomy care. You can depend on our experience, our product quality – and our quality of aftersales service. To order any of the products in this brochure or for details of other Kapitex products contact Kapitex on 01937 580211, or by fax on 01937 580796 or write to us at the address below.